

# Your Safety is our Priority

## European Standard EN 13795

### Part 1 - General Requirements

The first part describes the fundamental properties of surgical drapes, operating theatre clothing and clean air suits, and defines which requirements the product has to fulfil. Single-use and reusable products are subject to the same safety standards which apply for reusable products during the entire life cycle, thus also after every reprocessing. Products, which are intended for reuse, have to be provided with information on cleaning, disinfection, packing and methods of sterilisation and the number of reuses.

### Part 2 – Test Methods

The second part of the European Standard EN 13795 describes the test methods, under which the various properties of the material are tested. These tests stimulate everyday use of the product.

## Foliodrape® Disposable Surgical Drapes

Foliodrape® Disposable Surgical Drapes provide effective draping solutions for all types of surgery, exceeding EN 13795 standards and with a patented EAN traceability label. The choice of standard sets, single packed drapes and supplementary products ensures you can select the right product with suitable properties for every demand. Both Foliodrape® Comfort and Protect offer the highest level of fluid absorption across the whole drape. All drapes are lint and latex-free with anti-static properties. All draping sets are laid out in the order of use for quick and easy handling.



### Part 3 – Limiting Values

The third part defines the limiting values that the product has to achieve. As the risk of infection may vary depending on the duration of procedure and the volume of liquids, two compulsory performance levels are identified:

- **Standard Performance:** short duration and low volume of liquids
- **High Performance:** longer procedures and high volumes of liquids

In addition, the manufacturers have the possibility of defining critical (wound contact) or less critical (no wound contact).



helps healing.

**PAUL HARTMANN LTD**  
Heywood Distribution Park  
Pilsworth Road  
Heywood  
Lancashire  
OL10 2TT

Tel: 01706 363200  
Fax: 01706 363201

Web: [www.hartmann.co.uk](http://www.hartmann.co.uk)  
E-mail: [info@uk.hartmann.info](mailto:info@uk.hartmann.info)

## HARTMANN - Your partner in healthcare and hygiene

As a European leader in operating theatre products, HARTMANN offers a comprehensive range of disposable operating theatre drapes, gowns and accessories.

For further details, please request literature or visit the HARTMANN website at [www.hartmann.co.uk](http://www.hartmann.co.uk).



helps healing.

## Draping quality you can rely on

### Foliodrape® Comfort and Protect



## Foliodrape® Comfort for optimum comfort and safety

- In-patient operations are often time-consuming and involve large quantities of fluids.
- A surgical drape has to meet both the surgeon's and patient's requirements.
- Foliodrape® Comfort is designed not only to provide reliable protection against fluids but also to guarantee perfect fluid management and high skin comfort for the patient.



**Performance without compromise:**  
The proven Foliodrape® Comfort 3-layer drape.

### Optimum absorption strength:

The innovative Foliodrape® non-woven with 3D "orange skin structure" in the surface layer absorbs liquid quickly and reliably, is practically lint-free and offers good drapeability.

### Optimum barrier effect:

The pore-free, impermeable intermediate layer of polyethylene film is absolutely impermeable to fluids.

### Optimum comfort:


The bottom layer of non-woven(SMS) provides additional strength and stability, absorbs perspiration and rests with pleasant softness on the patient's skin.

### Benefits:

- Extra absorption capacity
- Reliable protection against bacteria
- Very drapable
- High patient comfort

## Foliodrape® Protect for an optimum combination of functionality and economy

- Short, surgical interventions in in-patient as well as out-patient settings put very specific demands on the surgical drape.
- Whilst infection control issues are paramount, choosing a product that meets the surgical requirements is economically appropriate.
- In addition to a unique 2-layer material, Foliodrape® Protect offers tailor-made surgical draping sets particularly suitable for minor surgical interventions.

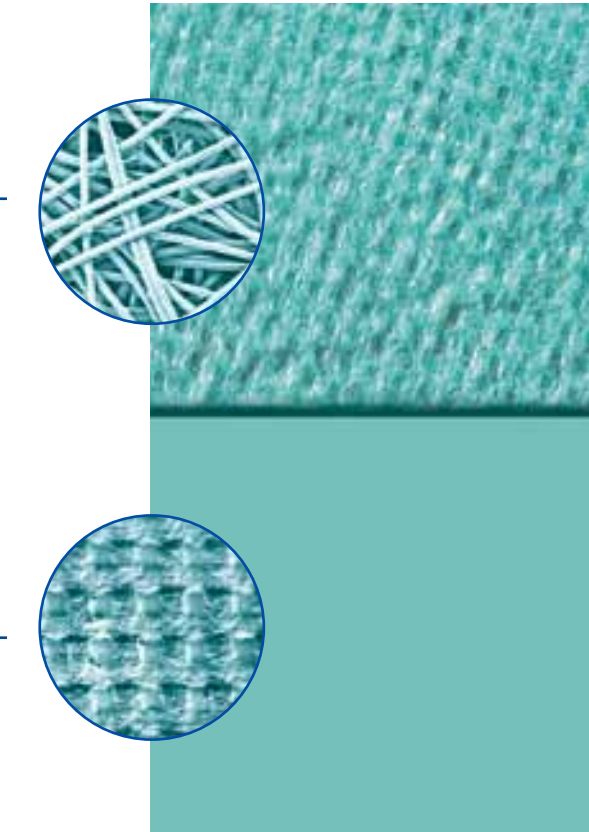

**Meets all demands on price and performance:**  
The Foliodrape® Protect 2-layer drape.

**Able to absorb any fluids:**  
The innovative Foliodrape® Protect fibre in the surface layer absorbs fluid quickly, is practically lint-free and drapes very well.

**Impermeable to fluids and bacteria:**  
Pore-free, impermeable polyethylene film

### Benefits:

- Excellent price-performance ratio
- Draping sets according to requirements
- Can be easily used by one person
- Very drapable