



HARTMANN



helps healing.

Zetuvit® Plus super absorbent wound dressing pad

For the treatment of superficial, heavily exuding wounds



Multi-layer structure

- Hydrophobic outer surface prevents adherence to the wound
- Hydrophilic cellulose fibres quickly pass exudate in to the absorbent core
- Green, water repellent, air permeable non woven layer protects against contamination

Super absorbent core

- Soft cellulose fluff blended with Super Absorbent Polymer (SAP) quickly absorbs and retains exudate
- Can be applied under compression therapy
- Soft consistency provides a pleasant cushioning effect



helps healing.

PAUL HARTMANN LTD
Heywood Distribution Park
Pilsworth Road
Heywood
Lancashire
OL10 2TT
Tel: 01706 363200
Fax: 01706 363201

Web: www.hartmann.co.uk
E-mail: info@uk.hartmann.info

Zetuvit® Plus super absorbent wound dressing pad

For the treatment of superficial, heavily exuding wounds

Zetuvit® Plus has a unique combination of absorption capacity, softness and cushioning. The blend of cellulose fluff and fluid retaining super absorbent particles (SAP) means that Zetuvit® Plus is particularly suitable for the treatment of superficial, heavily exuding, acute and chronic wounds.

In laboratory tests using protein containing test solution, Zetuvit® Plus exhibits the fastest absorption of exudate; the best absorption capacity under pressure and a free absorption capacity equal to that of the competition⁽¹⁾

Indications

Zetuvit® Plus is indicated for superficial, heavily exuding acute and chronic wounds.

Contraindications


None known.

⁽¹⁾Internal laboratory testing, data on file. February 2010.

Condition Dry



Condition Wet

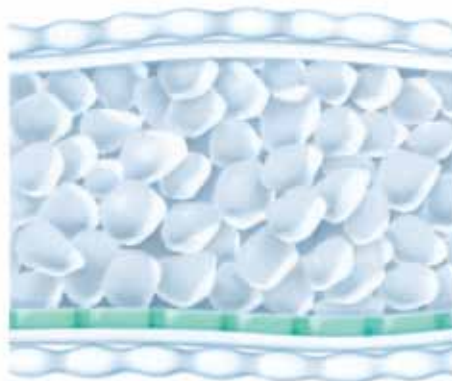

Zetuvit Plus provides excellent exudate management demonstrated by:

- + Rapid absorption of wound exudate**
Excess fluid is quickly absorbed and guided away into the core which considerably reduces the risk of exudate accumulation.
- + Excellent softness**
Zetuvit® Plus is extremely soft, has no sharp or protruding edges, providing additional relief for wounds and wound margins.
- + High-wearing comfort**
Due to the high proportion of fluff in the product, Zetuvit® Plus has an excellent padding effect right from the start.
- + Optimum cost-effectiveness**
Zetuvit® Plus is the economically sensible solution when treating heavily exuding wounds. Because of its excellent absorption capacity, fewer dressing changes and therefore less material are required.

Product	Size	Pack	Hartmann Code	NHS Code	PIP Code
Zetuvit® Plus	10cm x 10cm	Pack of 10	413710	EME046	356 7351
Super absorbent	10cm x 20cm	Pack of 10	413711	EME047	356 7385
wound dressing pad	15cm x 20cm	Pack of 10	413712	EME048	356 7369
	20cm x 25cm	Pack of 10	413713	EME049	356 7393
	20cm x 40cm	Pack of 5	413714	EME050	356 7377